



ICE SHAPING IV PRO
CRYOLIPOLYSIS SLIMMING

FREEZE YOUR FAT AWAY

One session. One hour.
No surgery. No Needles No downtime.




Treatment principle

The Cryolipolysis procedure safely delivers precisely controlled cooling to gently and effectively target the fat cells underneath the skin. Treated fat cells are crystalized (frozen), then die. Over time, your body naturally processes the fat and eliminates these dead cells, leaving a more sculpted you.



1

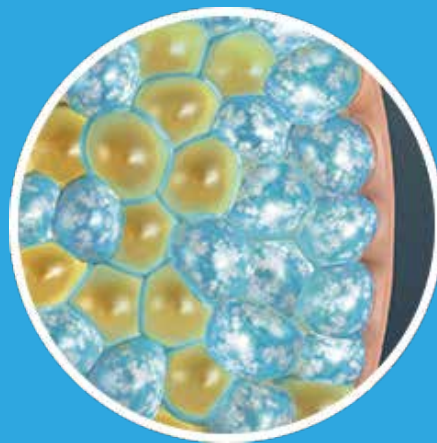

2



3




4

FAT REDUCTION VS WEIGHT LOSS

It may be confusing to hear that reducing fat is not the same as losing weight. It's all about what happens to fat CELLS. Losing weight makes fat cells smaller, losing weight does not make fat cells go away


The number of fat cells in our bodies typically becomes a fixed number sometime in our teenage years. Losing or gaining weight generally does not increase or decrease the number of fat cells. Instead, it changes the size of fat cells.

Fat Cells After Typical Weight Loss

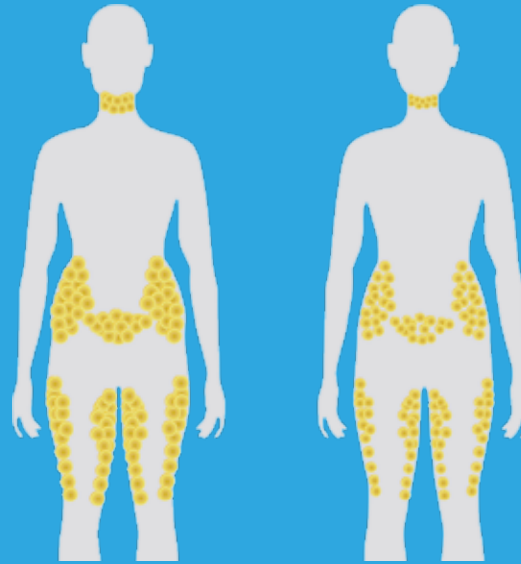

After dieting and weight loss, fat cells get smaller but the number of fat cells stays the same. When weight is gained, fat cells get bigger.

Fat Cells After The Cryolipolysis Procedure

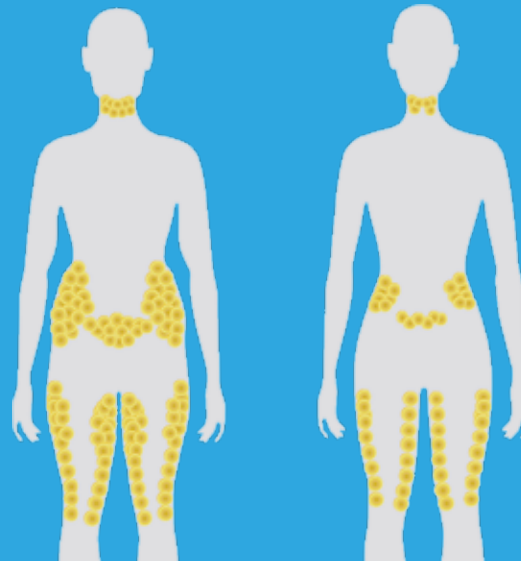

The Cryolipolysis fat-freezing procedure actually reduces the number of fat cells in treated areas by about 20% to 25%. These fat cells are naturally processed and eliminated, so they won't migrate to another area of your body.






Advantages

- Four different sizes handpieces optional, easily treatment for whole body.
- Four handpieces working simultaneous, save staff cost, save treatment time.
- Unique double chin treatment handpiece, lead the latest technology in the world.




- Super strong cooling system ensure the minimum treatment temperature reduce to -15°C
- Non-surgery, no risk, easily intelligentize and humanization operation
- No need recovery, no downtime, come back to work immediate.

The IEC SHAPING IV PRO has 4 different sizes handpieces optional, with these 4 handpieces, you can do nearly whole body fat freezing.

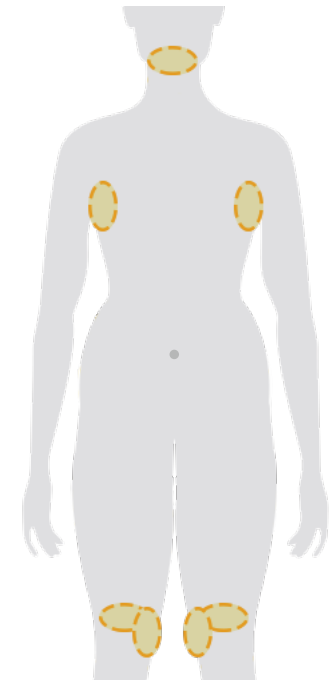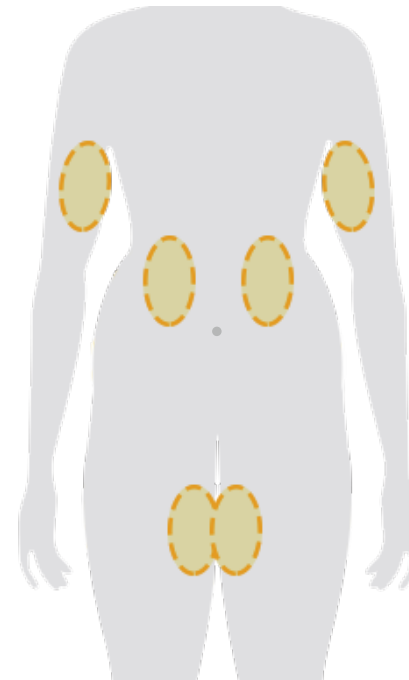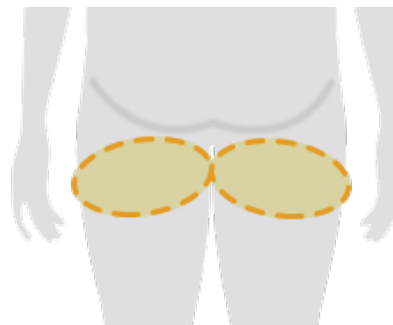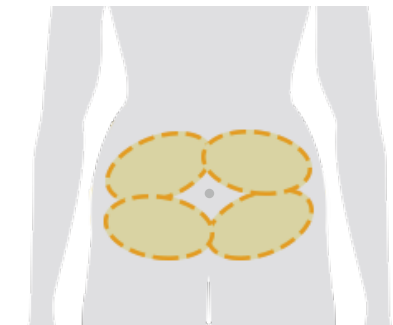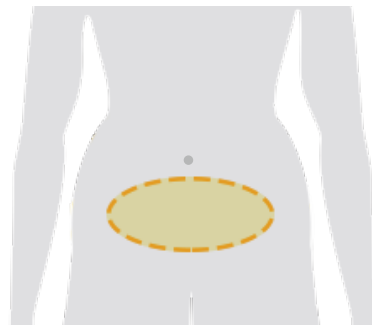



The CryoMini can not only do double chin treatment, but also do all the small areas fat freezing.

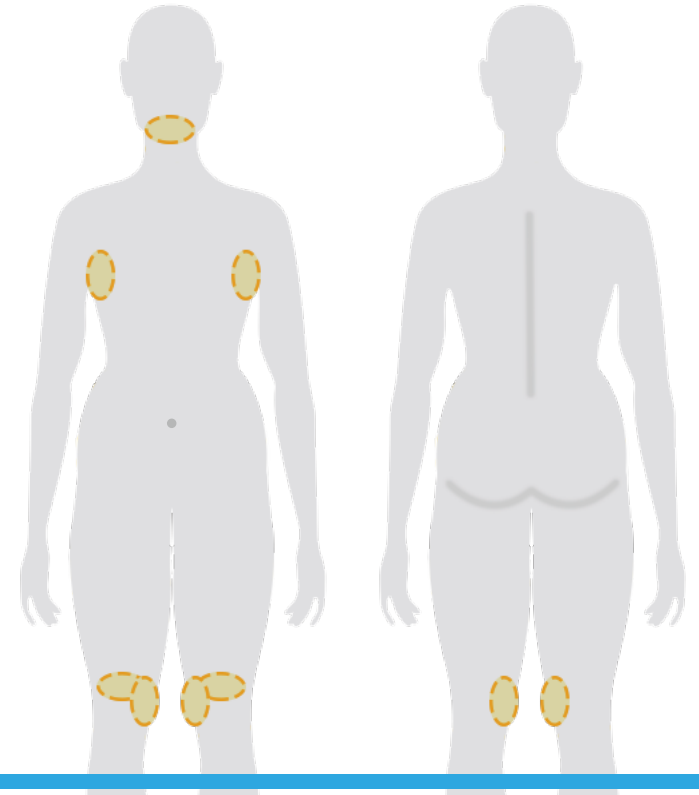

**SAY GOODBYE TO YOUR
DOUBLE CHIN**



Specification



Handle number	4 handles available
Handle screen	5 inch touch screen
Vacuum	0kpa-100kpa
Handle size	Big handle A: 220 * 76*125mm Middle handle B1: 160 * 56*65mm Middle handle B2: 160 * 56*65mm Small handle C: 125 * 45*70mm
Screen	10.4 inch touch screen
Cooling temperature	5 --15℃
Cooling system	Wind cooling + water cooling+ semiconductor cooling
input power	2500W
Output power	1600W
Dimension	1160*508*620mm
Net weight	45kg
Fuse specification	Ø5×25 15A
Voltage	AC220V±10%,10A;50HZ/ AC110V±10%10A,60HZ

Handle Advantages

- Suction cup is made of soft medical use silica gel; the whole treatment process is comfortable



- Eight Peltier cooling plates are used inside the handle; the lowest treatment temperature can reach -15°C



- 5.0 inch touch screen, clearly displaying the treatment parameters.
- Three minutes heating first makes a better blood circulation of the treated area to achieve the best treatment results.
- Two temperature sensor in each side, ensure the safety of temperature control.



Before & After



Different Angles

ICE SHAPING IV PRO

